



Illumination Equipment

Light Sources – Xenon Illuminators.....	M-1
Light Source Replacement Parts	M-2
Fiber Optic Light Cables and Adaptors	M-3
Compatibility Matrix	M-5

Light Sources – Xenon Illuminators

Smith & Nephew DYONICS® Light Source provides bright white light for natural color reproduction. Features include a multi-port light cable receptacle for use with fiber optic light cables with DYONICS, Olympus®, Storz®, or ACMI® source ends.

300XL Xenon Light Source, 300W



REF	Description
7207407*	300XL System, English (U.S. power cord), 100-240 VAC
7207399	300XL System, Germany, 100-240 VAC
7207401	300XL System, France, 100-240 VAC
7207402	300XL System, Italy, 100-240 VAC
7207400	300XL System, Spain, 100-240 VAC
7207403	300XL System, Swedish, 100-240 VAC
7207404	300XL System, Dutch, 100-240 VAC
7207405	300XL System, English (U.K. power cord), 100-240 VAC
7207406	300XL System, English (Continental power cord), 100-240 VAC

CONDOR®-Enabled 300XL Xenon Light Source, 300W




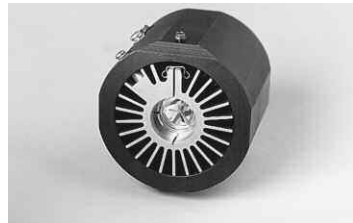

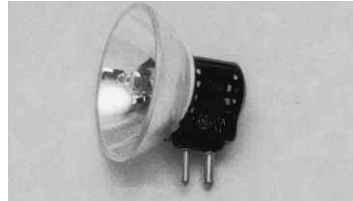
Adds CONDOR compatibility to the Smith & Nephew 300XL light sources.

7210261	300XL CONDOR-Enabled Light Source, English (U.S. power cord)
7210263	300XL CONDOR-Enabled Light Source, English (U.K. power cord)
7210265	300XL CONDOR-Enabled Light Source, International System

* Item shown

Light Source Replacement Parts

Replacement Lamps


REF	Description
7207398*	Replacement Xenon Lamp Module, 300W, for 300XL Xenon Light Source, 500-hour minimum life
	
72200529*	Replacement Xenon Lamp Module, 180W, for 420XL Light Source (72200510)
	
7205351	Replacement Xenon XL Lamp Module, 180W, for Xenon Light Source, 500-hour minimum life
	
7204741*	Replacement Xenon Lamp Module, 100W, for Xenon Light Source, 500-hour minimum life
	
4248*	Replacement Xenon Lamp Assembly for DYOBRITE® 3000 Light Source (REF 4243)
	
2974*	Replacement Metal Halide Arc Lamp with bracket, 250W for AutoBrite, MicroBrite, or ABI II II Illuminators, 250-hour life
4550	Replacement Arc Lamp High-Voltage Cable Assembly for AutoBrite Illuminator II
	
7204563*	Replacement Lamp for DYOSITE® II, 150 W, 100-hour life
7209213	Replacement Lamp for DV 111, 60 W, 500-hour life

* Item shown

Fiber Optic Light Cables and Adaptors

Gemini® Universal 4 mm Fiber Optic Light Cables without Adaptor

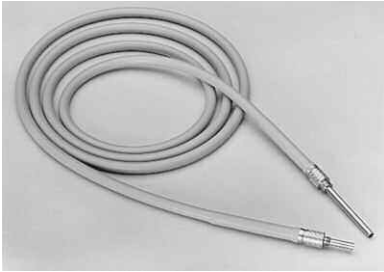For use with Smith & Nephew 4 mm and 2.7 mm arthroscopes.



REF	Description
2985*	4.0 mm x 8 ft
7205180	4.0 mm x 10 ft
7205179	4.0 mm x 12 ft

Use only REF 2146 and REF 2147 Versitip Adaptors for connection to Smith & Nephew Light Sources and endoscopes.


Gemini® Universal 5 mm Fiber Optic Light Cables without Adaptor

For use with all arthroscopes and laparoscopes.

REF	Description
2140	5.0 mm x 8 ft
7205178	5.0 mm x 10 ft
7205177	5.0 mm x 12 ft

Use only REF 2146 and REF 2147 Versitip Adaptors for connection to Smith & Nephew Light Sources and endoscopes.

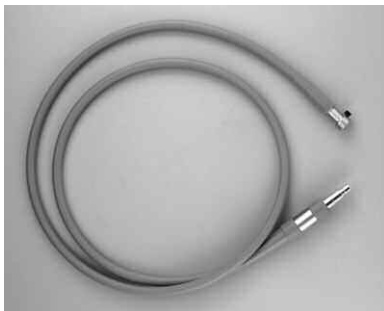Gemini® Universal Fiber Optic Light Cables with Adaptors



REF	Description
7208329	6.0 mm x 10 ft, with permanently fixed Smith & Nephew/Wolf® source and instrument end adaptors. For use with all 10 mm laparoscopes. Not for use with arthroscopes or 5.0 mm laparoscopes.






3004	2.0 mm x 8 ft, with Smith & Nephew/Wolf adaptors. For use with 1.9 mm arthroscopes (with REF 2147 in package).
------	--

* Item shown

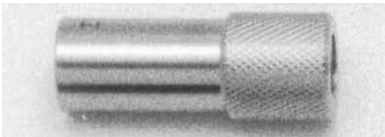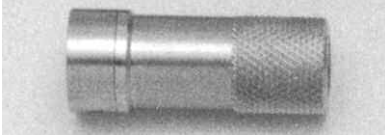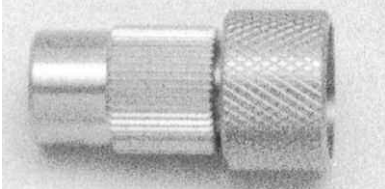Fiber Optic Light Cables and Adaptors




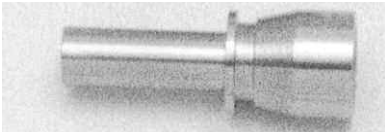
Versitip Adaptors – Instrument End

For Smith & Nephew Gemini® Universal Fiber Optic Light Cables.

	REF	Description
	2141	ACMI® Instrument ("AI") Versitip Adaptor for 4.0 mm and 5.0 mm light guides Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cables with ACMI endoscopes.
	2146	Smith & Nephew/Wolf® Instrument ("WI") Versitip Adaptor Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cables with Wolf and Smith & Nephew endoscopes.
	2143	Storz®/Olympus® Instrument ("SI") Versitip Adaptor Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cables with Storz or Olympus endoscopes.

Versitip Adaptors – Light Source End

For Smith & Nephew Gemini Universal Fiber Optic Light Cables.

	2142	ACMI Source ("AS") Versitip Adaptor Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cables with ACMI Light Sources. Not for use with Smith & Nephew light sources.
	2147	Smith & Nephew/Wolf Source ("WS") Versitip Adaptor Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cable with Smith & Nephew or Wolf light sources.
	2145	Olympus Source ("OS") Versitip Adaptor Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cable with Olympus light sources. Not for use with Smith & Nephew light sources.
	2144	Storz Source ("SS") Versitip Adaptor Intended to interface Smith & Nephew Gemini Universal Fiber Optic Light Cable with Storz light sources. For use with Smith & Nephew 420XL Light Source.

Compatibility Matrix

To connect a Smith & Nephew light guide to a Smith & Nephew light source:

Smith & Nephew Light Guide	Smith & Nephew Versitip Adaptor Light Source End	Light Source Vendor	Connect to Light Source Input Port labeled:
7208329	NONE	Smith & Nephew	Wolf®/Smith & Nephew
3004	2147	↓	↓
2985	2147		
7205180			
7205178			
7205179			
7205177			

To connect a Smith & Nephew light guide to another vendor's light source:

Smith & Nephew Light Guide	Smith & Nephew Versitip Adaptor Light Source End	Light Source Vendor	
2985	2142	ACMI®	
2140	2147	Wolf	
7205180	2145	Olympus®	
7205178	2144	Storz®	
7205179	↓	↓	
7205177			

To connect another vendor's light guide to a Smith & Nephew light source:

Light Guide Vendor	Smith & Nephew Versitip Adaptor Light Source End	Light Source Vendor	Connect to Light Source Input Port labeled:
ACMI	NONE	Smith & Nephew	ACMI
Wolf	Insert in corresponding receptacle on front panel labeled with vendor's name,	↓	Wolf/Smith & Nephew
Olympus			Olympus
Storz			Storz

